



Why Choose Klercide 70/30 IPA Sterile Spray 1L - Case of 6

Sold in cases of 6

Sterile, DI, ready to use, 1 litre (WFI version available)

70% isopropyl alcohol (IPA)

Double bagged & 0.2 micron filtered

SteriShield delivery system

Formerly - Premier Klercide 70/30 Sterile IPA Spray 1L

Very suitable for a number of applications in GMP Grades A/B/C/D

Ecolab code: 6600946 formally 3078480

A dangerous goods shipping surcharge will be applied to this product for shipments leaving the UK mainland

More Info

Klercide 70|30 IPA is a blend of 70% v/v Isopropyl Alcohol with deionized water.


The alcohol blend is 0.2 micron filtered, filled and double bagged in a Grade C (ISO class 7) cleanroom, before irradiation using a validated process at no less than 25 kGy to ensure a spore free alcohol.

All trigger spray bottles incorporate the patented SteriShield Delivery System and an adjustable trigger spray. The unique way in which the SteriShield Delivery System works as a closed system, ensures that the sterility of the contents is preserved throughout use.

Advantages

Easy to use - convenient range of comfortable and easy to use trigger sprays is available, with a lightweight option of 500ml. The large droplet size reduces the risk of inhalation.

Cost effective - all the sterile alcohol can be dispensed from the bottle, so there is no wastage. Validation of the SDS demonstrates that the contents remain sterile. An in-use shelf life of three months is supported with the data, enabling

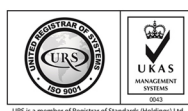

www.cleanroomshop.com

Tel: + 44 (0) 1524 813022

info@cleanroomshop.com

Fax: + 44 (0) 1524 811589

enabling thorough wetting of the surface so that disinfection procedures are effective.

More environmentally friendly - the product contains no propellants so does not require any special disposal procedures.

Product ordering

Format: 1 Litre SDS Trigger Spray

Case quantity: 6 x 1L

Ecolab code: 6600946 formally 3078480

[Safety datasheet \(PDF\)](#)

Please note: A dangerous goods surcharge may apply - more information can be found on our delivery information page.

Product Icon Attributes

Suitable for GMP A/B/C/D, Sterile, DI, 70/30 IPA

Brand

Ecolab

Type

SteriShield delivery system

Lead time (if out of stock)

10-15 days

Pack Size

Case of 6

Materials

70% isopropyl alcohol

Size

1L

Colour

Colourless

Odour

Alcohol-like

Clarity

Clear

Water grade

DI

Applications

Very suitable for a number of applications in Grades A/B/C/D areas - please enquire for full details



www.cleanroomshop.com
Tel: + 44 (0) 1524 813022

info@cleanroomshop.com
Fax: + 44 (0) 1524 811589

Production

Klercide 70 | 30 IPA blended with deionized water is 0.2 micron filtered then filled and double bagged in a Grade C cleanroom

Packaging

Double bagged

Shelf life

2 years from date of manufacture

Dilution

Ready to use

Cleanroom class

GMP Grade A, GMP Grade B, GMP Grade C, GMP Grade D

Handling & storage

Store in cool well ventilated areas. Keep container tightly closed. Must only be kept in original packaging.

Instructions for use

Ready to use so can be applied directly onto the surface to be sanitised. Always ensure the whole surface is in contact with the liquid. For optimum results wipe the surface to physically remove any contaminants.

Country of origin

United Kingdom

Sterility

Sterile

Efficacy

CEN Test Data available on request. Certificates of analysis, irradiation, quantitative endotoxin levels < 0.25EU/ml and sterility available

MSDS


Available to download on the Cleanroomshop product page

PPE Requirements

Ensure that there is sufficient ventilation to the area. Nitrile, PVC or Neoprene gloves for hand protection. Eye glasses for eye protection and protective clothing for overall protection.

Dangerous Goods

Yes

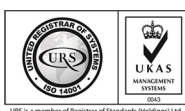

www.cleanroomshop.com
Tel: + 44 (0) 1524 813022

info@cleanroomshop.com
Fax: + 44 (0) 1524 811589