



Why Choose Klercide 70/30 IPA Sterile Capped 4 x 5L

4 x 5L bottles per case, sold by the case

Sterile, DI, ready to use, 5 litre

70% isopropyl alcohol (IPA)

Double bagged & 0.2 micron filtered

SteriShield delivery system

Great value - to be decanted

Very suitable for a number of applications in Grades A/B/C/D areas

Formerly - Premier Klercide 70/30 Sterile IPA Capped 5L

Ecolab code: 3078460

A dangerous goods shipping surcharge will be applied to this product for shipments leaving the UK mainland

More Info

Klercide 70|30 IPA is a blend of 70% v/v Isopropyl Alcohol with deionized water.

The alcohol blend is 0.2 micron filtered, filled and double bagged in a Grade C (ISO class 7) cleanroom, before irradiation using a validated process at no less than 25 kGy to ensure a spore free alcohol.

Product ordering

Format: 5 litre capped

Case quantity: 4 x 5L

Ecolab code: 3078460

[Safety datasheet](#) (PDF)

Please note: A dangerous goods surcharge may apply - more information can be found on our delivery information page.

Product Icon Attributes

Suitable for GMP A/B/C/D, Sterile, DI, 70/30 IPA

Brand

Ecolab

Type

Capped Bottle

Lead time (if out of stock)

10-15 days

Pack Size

4 x 5L bottles per case, sold by the case

Materials

70% isopropyl alcohol

Size

5L

Colour

Colourless

Odour

Alcohol-like

Clarity

Clear

Water grade

DI

Applications

Very suitable for a number of applications in Grades A/B/C/D areas - please enquire for full details

Production

Klercide 70 | 30 IPA blended with deionized water is 0.2 micron filtered then filled and double bagged in a Grade C cleanroom

Country of origin

United Kingdom

Packaging

Double bagged

Sterility

Sterile

Shelf life

2 years from date of manufacture

Efficacy

CEN Test Data available on request. Certificates of analysis, irradiation, quantitative endotoxin levels < 0.25EU/ml and sterility available

Dilution

Ready to use

MSDS

Available to download on the Cleanroomshop product page

Cleanroom class

GMP Grade A, GMP Grade B, GMP Grade C, GMP Grade D

PPE Requirements

Ensure that there is sufficient ventilation to the area. Nitrile, PVC or Neoprene gloves for hand protection. Eye glasses for eye protection and protective clothing for overall protection.

Handling & storage


Store in cool well ventilated areas. Keep container tightly closed. Must only be kept in original packaging.

Dangerous Goods

Yes

Instructions for use

Ready to use so can be applied directly onto the surface to be sanitised. Always ensure the whole surface is in contact with the liquid. For optimum results wipe the surface to physically remove any contaminants.

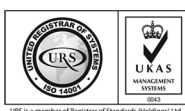

www.cleanroomshop.com
Tel: + 44 (0) 1524 813022

info@cleanroomshop.com
Fax: + 44 (0) 1524 811589